CHORUS

INTRINSIC FACTOR

REF 86090

REF 86090/12

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INSTRUCTIONS FOR USE

CHORUS INTRINSIC FACTOR

For the semiquantitative determination of IgG antibodies anti-Intrinsic Factor

For *In Vitro* Diagnostic Use Only

1. INTENDED USE

Immunoenzymatic method for the semiquantitative determination of IgG class antibodies against Intrinsic Factor in human serum, using a disposable device applied on the chorus and Chorus TRIO instruments.

2. INTRODUCTION

Pernicious Anaemia is the end stage of an autoimmune gastritis (type A gastritis) resulting in destruction of the gastric mucosa. A small proportion (10-15%) of patients with autoimmune gastritis develop Pernicious Anaemia. Pernicious Anaemia is the most common cause of vitamin B12 deficiency in Western populations. It is caused by a lack of Intrinsic Factor, a glycoprotein needed to absorb vitamin B12 from the gastrointestinal tract. Intrinsic factor is produced by gastric parietal cells. Vitamin B12, in turn, is necessary for the formation of red blood cells.

Pernicious Anaemia usually does not appear before the age of 30. The average age at diagnosis is 60 years. In fact, one recent study revealed that nearly 2% of individuals over 60 years old suffer from Pernicious Anaemia. Furthermore, slightly more women than men are affected. The disease can affect all ethnic groups, but occurs more often among people of Scandinavian and Northern European descent.

Autoantibodies in Pernicious Anaemia target parietal cells and Intrinsic Factor. Antibodies against parietal cells are 80-90% sensitive but can be detected also in up to 5% of the healthy population. Antibodies against Intrinsic Factor show a sensitivity of 50-70% with a specificity of 100% in a population of healthy blood donors.

3. PRINCIPLE OF THE METHOD

The Chorus Intrinsic Factor device is ready to use for the detection of IgG antibodies against Intrinsic Factor, in the Chorus/Chorus TRIO instruments.

The test is based on the ELISA principle (Enzyme linked Immunosorbent Assay). The antigen is bound to the solid phase. The specific immunoglobulins are bound to the antigen through incubation with diluted human serum. After washing to eliminate the proteins which have not reacted, incubation is

performed with the conjugate, composed of anti-human immunoglobulins antibodies conjugated to horse radish peroxidase.

The unbound conjugate is eliminated and the peroxidase substrate is added. The colour which develops is proportional to the concentration of specific antibodies present in the serum sample.

The disposable devices contain all the reagents to perform the test in the Chorus/Chorus TRIO instruments.

The results are expressed in Arbitrary Units (AU/ml).

4. WARNINGS AND PRECAUTIONS

FOR *IN VITRO* DIAGNOSTIC USE ONLY

This kit contains materials of human origin which have been tested and gave a negative response by FDA-approved methods for the presence of HBsAg and for anti-HIV-1, anti-HIV-2 and anti-HCV antibodies. As no diagnostic test can offer a complete guarantee regarding the absence of infective agents, all material of human origin must be handled as potentially infectious. All precautions normally adopted in laboratory practice should be followed when handling reagents and samples.

Waste disposal: serum samples, calibrators and strips, once used, must be treated as infectious residuals and eliminated according to law.

Health and Safety Information

1. Do not pipette by mouth.
2. Wear disposable gloves and eye protection while handling specimens.
3. Wash hands thoroughly after placing the devices in the Chorus/Chorus TRIO instrument.
4. Consult the relative Material Safety Data Sheet (available on request) for all the information on safety concerning the reagents contained in the kit.
5. Neutralized acids and other liquid waste should be decontaminated by adding a sufficient volume of sodium hypochlorite to obtain a final concentration of at least 1%. A 30 minute exposure to 1% sodium hypochlorite may be necessary to ensure effective decontamination.
6. Spillage of potentially infectious materials should be removed immediately with adsorbent paper tissue and the contaminated area swabbed with, for example, 1% sodium hypochlorite before work is continued. Sodium hypochlorite should not be used on acid-containing spills unless the spill area is first wiped dry. Materials used to clean spills, including gloves, should be disposed of as potentially biohazardous waste. Do not autoclave materials containing sodium hypochlorite.

Analytical Precautions

Bring the devices to room temperature (18-30°C) before use; use within 60 min.

1. **Discard devices which show the substrate (well 4) blue colored.**
2. Adding the sample into the well verify that it is perfectly distributed on the bottom.
3. Check for the presence of the reagents in the device and that the device is not damaged. Do not use devices which are lacking a reagent and/or present foreign bodies in the reaction well when visually inspected.
4. The devices are for use with the Chorus/Chorus TRIO instrument; the Instructions for Use and the Instrument Operating Manual must be carefully followed.
5. Check that the Chorus/Chorus TRIO instrument is set up correctly (see Operating Manual).
6. Do not alter the bar code placed on the handle of the device in order to allow correct reading by the instrument.
7. Avoid using self-defrosting freezers for the storage of the samples.
8. Defective barcodes can be inserted manually in the instrument (see Operating Manual).
9. Do not expose the devices to strong light or to hypochlorite vapors during storage and use.
10. The use of strongly hemolyzed, lipemic, icteric samples, of serum not completely coagulated or of samples presenting microbial contamination may all constitute sources of error.
11. Do not use the device after the expiry date.
12. **Make sure that the instrument is connected to the Washing Buffer Autoimmunity (Ref. 86004).**

5. KIT COMPOSITION AND REAGENT PREPARATION

The kit is sufficient for 36 tests (REF 86090).


The kit is sufficient for 12 tests (REF 86090/12).

DD DEVICES

6 packages each containing 6 devices (REF 86090).

2 packages each containing 6 devices (REF 86090/12).

Description:

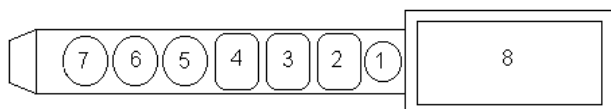Position 8: Space for application of bar code label

Position 7: Empty

Position 6: MICROPLATE WELL

Coated with recombinant Intrinsic Factor

Position 5: Uncoated MICROPLATE WELL

Position 4: TMB SUBSTRATE

Contents: Tetramethylbenzidine 0.26 mg/mL and H₂O₂ 0.01% stabilized in 0.05 mol/L citrate buffer (pH 3.8)

Position 3: SAMPLE DILUENT

Contents: Proteic solution in Tris, saline, with sodium azide

Positon 2: CONJUGATE

Contents: anti-human IgG monoclonal antibodies labeled with horse radish peroxidase.

Position 1: EMPTY WELL

in which undiluted serum must be added

Use: equilibrate a package at room temperature, open the package and remove the required devices; replace the others in the bag with the silica gel, expel the air and **seal** by pressing the closure. Store at 2-8°C.

CALIBRATOR CALIBRATOR 1 x 0.175 mL

Contents: Diluted human serum. Liquid, ready for use.

CONTROL + POSITIVE CONTROL 1 x 0.425 mL

Contents: Diluted human serum. Liquid, ready for use.

MATERIALS REQUIRED BUT NOT PROVIDED:

- WASHING BUFFER AUTOIMMUNITY (REF 86004)
- CLEANING SOLUTION 2000 (REF 83609)
- SANITIZING SOLUTION (REF 83604 - 83608)
- CHORUS NEGATIVE CONTROL/SAMPLE DILUENT (REF 83607)
- Chorus/Chorus TRIO Instrument
- Distilled or deionised water
- Normal laboratory glassware: cylinders, test-tubes etc.
- Micropipettes for the accurate collection of 50-200 µl solution
- Disposable gloves
- Sodium Hypochlorite solution (5%)
- Containers for collection of potentially infectious materials

6. STORAGE AND STABILITY OF REAGENTS

Reagents must be stored at 2/8°C. In the case of storage at an incorrect temperature the calibration must be repeated and the run validated using the control serum (see section 9, Test validation).

The expiry date is printed on each component and on the kit label.

Reagents have a limited stability after opening:

DEVICES	8 weeks at 2/8°C
CALIBRATOR	8 weeks at 2/8°C
POSITIVE CONTROL	8 weeks at 2/8°C

7. SPECIMEN COLLECTION AND STORAGE

The sample is composed of serum collected in the normal manner from the vein and handled with all precautions dictated by good laboratory practice.

Possible consequences, in case of use of other biological liquids, are not known.

The fresh serum may be stored for 4 days at 2/8°C, or frozen for longer periods at -20°C, and can be thawed a maximum of 3 times.

Do not keep the samples in auto-defrosting freezers. Defrosted samples must be shaken carefully before use. Heat-inactivation can rise to erroneous results. The quality of the sample can be seriously affected by microbial contamination which leads to erroneous results.

8. ASSAY PROCEDURE

1. Open the package (on the side containing the pressure-closure), remove the number of devices required and seal the rest in the bag after expelling the air.
2. Check the state of the device according to the indications reported in chapter 4, Analytical Precautions.
3. Dispense 50 µl of undiluted test serum in well no. 1 of each device; at each change of batch, use a device for the calibrator.
4. Place the devices in the Chorus/Chorus TRIO instrument. Perform the calibration (if necessary) and the test as reported in the instrument Operating Manual.

9. TEST VALIDATION

Use the control serum to check the validity of the results obtained. It should be used as reported in the instrument Operating Manual. If the instrument signals that the control serum has a value outside the acceptable range, the calibration must be repeated. The previous results will be automatically corrected.

If the result of the control serum continues to be outside the acceptable range, contact the Scientific Support.

Tel: 0039 0577 319554
Fax: 0039 0577 366605
email: scientificsupport@diesse.it

10. INTERPRETATION OF THE RESULTS

The Chorus/Chorus TRIO instrument expresses the result in Arbitrary Units (AU/ml) calculated on the basis of a lot-dependent graph stored in the instrument.

The test on the examined serum can be interpreted as follows:

POSITIVE: when the result is > 18 AU/ml

NEGATIVE: when the result is < 12 AU/ml

DOUBTFUL/EQUIVOCAL: for all values between 12 and 18 AU/ml

If the result is doubtful/equivocal, repeat the test. If it remains doubtful/equivocal, collect a new serum sample.

11. LIMITATIONS

All the values obtained require a careful interpretation that must consider other indicators relative to the patient.

The test, indeed, can not be used alone for a clinical diagnosis and the test result should be evaluated together with the patient history and other clinical diagnostic evaluation.

12. CALIBRATION RANGE

Calibration range 3.0-300.0 AU/ml.

For samples > 300.0 AU/ml retest the diluted sample in the Negative Control/Sample Diluent (PF83607-not supplied with the kit).

13. ANALYTICAL SENSITIVITY

The analytical sensitivity of this kit has been found at 3.0 AU/ml.

14. SENSITIVITY AND SPECIFICITY

The microplates are coated with recombinant human gastric intrinsic factor. No crossreactivities to other autoantigens have been found. Antibodies against intrinsic factor show a diagnostic sensitivity of 50-70% for Pernicious Anaemia and show a specificity of 100% in a population of healthy blood donors.

15. LINEARITY

Chosen sera have been tested with this kit and found to dilute linearly. However, due to the heterogeneous nature of human autoantibodies there might be samples that do not follow this rule.

Sample N°.	Dilution factor	Measured concentration (AU/ml)	Expected concentration (AU/ml)	Recovery (%)
1	Undiluted	160	163	98.1
	1/2	79	81	97.5
	1/4	37	41	90.2
	1/8	19	20	95
2	Undiluted	101	104	97.1
	1/2	53	52	101.9
	1/4	28	26	107.7
	1/8	12	13	92.3

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DIESSE Diagnostica Senese S.p.A.